

⚠ Safety precautions

Important note:

The terahertz energy machine acts on the human body through high-intensity electromagnetic energy field, regulates health, and controls high-intensity energy must cooperate with scientific and accurate operation.

1. Please read the "Safety Precautions" standard operation, to avoid improper operation of electric shock, burning and other adverse reactions, wish every user can get a safe and good product experience!

2. The three safety levels: [Danger] [Warning] [Attention] are very important, must comply with!

⚠ Danger:



Make sure your feet closely contact with the energy board



Don't touch your feet or the energy board when using the machine



Forbidden to use it with other device which using electromagnetic technology



No contact for both feet (including knees) when using

⚠ Warning:



If you feel uncomfortable during use, please decrease intensity or stop using



It's forbidden to place anything on the energy board



Please keep the machine away from children, use it by people with ability



Please put both feet on the energy board, don't leave one of your foot suddenly from the energy board when using

! Warning:



In case of power failure during use, etc. When the situation, should be timely disconnect the power supply



Please do not use any metal objects to touch the energy board



Do not use under high temperature or high humidity environment



Please take off all metal objects before use

Usage taboo

1. Individuals with low blood sugar. Patients with diabetes and malignant tumors.

2. Those with high blood pressure (systolic pressure ≥ 165 or diastolic pressure ≥ 110) or symptoms of cerebral hemorrhage.

3. Individuals with high fever, hemorrhagic diseases, bleeding tendencies, trends of critical organ failure, or symptoms of cardiopulmonary dysfunction.

4. Patients with heart disease and those with electronic devices implanted in the body such as cardiac pacemakers, as well as those with metal objects like steel plates or orthopedic braces implanted in the body.

5. Patients with heart and kidney failure, severe cardiac, liver, lung, or kidney failure, those who are extremely frail, and patients with serious illnesses.

6. Patients who have had a heart attack, stroke, or cerebral infarction and those who have had a stroke less than six months ago.

7. Pregnant women, women in the first year of breastfeeding, women during menstruation, infants and young children, patients with epilepsy, mental illness, and individuals under 18 years of age are prohibited from use; those over 60 years of age should use it with assistance.

8. Individuals who have undergone major surgery within the past three months, those with inflamed surgical wounds, and those with bleeding or a tendency to bleed.

9. People who are insensitive to temperature, alcohol consumers, and those unable to care for themselves independently.

Safety Guide

Thank you for choosing our product! To prevent accidents, please read the user manual thoroughly before using this product. If you are transferring the product to someone else, ensure the manual is included.

1. Before use, check that the voltage matches the one indicated on the device.

2. Familiarize yourself with the operating methods and precautions before use.

3. Do not spill any liquids on the device or allow them to enter the device to avoid the risk of electric shock. Turn off and unplug the device before cleaning to prevent electric shock.

4. Do not use this device in conjunction with other devices during therapy, and do not stand on the device.

5. Use only the provided or recommended accessories. Using other types of accessories may damage the device and affect its performance and safety.

6. Users should limit each session to 30 minutes, with a maximum of 40 minutes per use. Ensure both feet are placed on the pedals in the standard position before starting use. Allow a pause of more than 5 minutes in the off state between two consecutive uses of the device. Do not keep your knees together, and avoid touching the pedals and metal panels with your hands during use to prevent burns.

7. Do not apply external impacts to the machine, and do not modify, disassemble, or repair the device, as this will void the warranty.


8. If your knees are exposed, do not keep them together to avoid burns. It is normal for athletes' foot patients to experience blistering, indicating a good reaction. However, there is generally no need for concern about injury.

9. During use, both feet must be in close contact with the pedals, avoiding lifting one foot, and ensure balanced lifting to prevent concentrated current or rapid heating at a single point, which may cause momentary high heat or a tingling sensation.

Disclaimer: This device is not a medical instrument and is strictly prohibited for high-risk individuals! Our company will not be responsible for any personal or property damage caused by improper use or storage by the customer.

USER MANUAL


TERAHERTZ ENERGY MACHINE

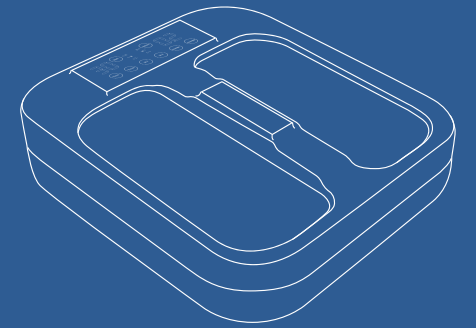


MODEL NO: GY-SQT-02

Description of components


Product position function

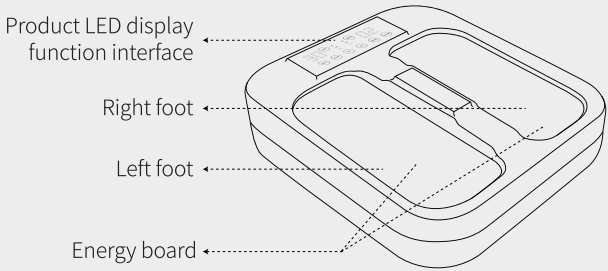


Product LED display function interface

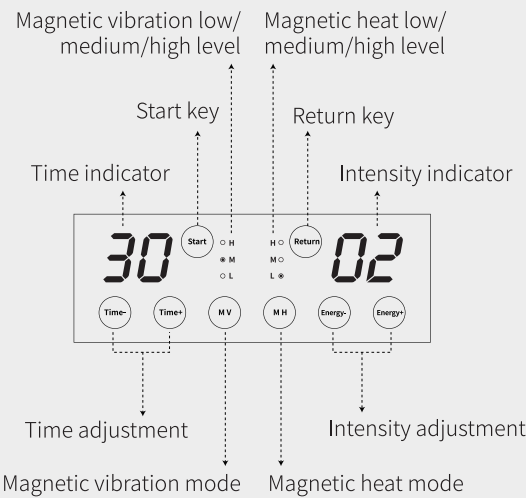

Product remote control function interface

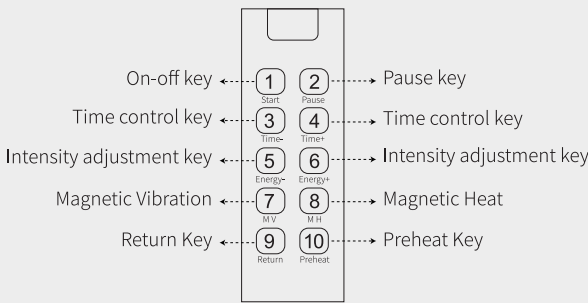

Product back position function

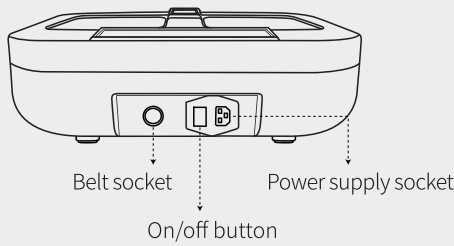

Pairing the remote controlwith the main unit:


When the remote control and the host are not responding, they need to be paired.

First, press and hold the "Learn" button on the remote control, and then turn on the power switch of the host. Once you hear a beep, beep, beep. Release the button and the device is ready to use.

Tip:

This device must have both feet placed on the energy board to start working. If only one foot is placed, it will not work and an alarm will sound. The device will also not operate if both feet are not properly positioned on the energy board.

Product Parameters

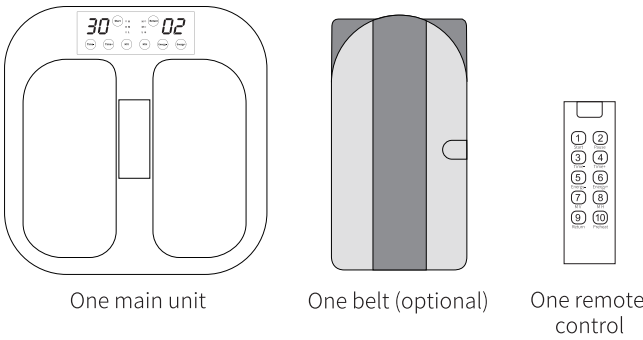

Product Name	Terahertz Energy Machine	Model	GY-SQT-02
Voltage	AC 110V <input type="checkbox"/> AC 220V <input type="checkbox"/>	Frequency	50Hz <input type="checkbox"/> 60Hz <input type="checkbox"/>
Product Size	41x38x10cm/16x15x3.5in	Power	300W
Package Size	50x42x19cm/20x16.5x7in	N.W	3.15KG
Packing List	Main unit, Belt (optional), Remote controller, User manual, power line, Two spare fuses		

Special caution: Patients with severe hypertension or heart conditions should be aware that using this device may increase blood circulation. It must be used under supervision, and the output intensity and duration should be adjusted gradually based on the user's physical response.

For first-time users, it's recommended not to exceed 30 minutes of using, the energy intensity should be set started Low. After using the device for 7 days, can increase the duration and energy intensity according to the user's using feeling.

Instructions for use

1. Connect to the power source.
2. Activate Bluetooth on your mobile phone, search for "Rotary Magnetic Wave Instrument," and click to connect the device via Bluetooth. This will allow the device to play audio content from your phone.
3. Position both feet on the energy board according to the standard posture; knees should not be touching. Avoid touching the energy board with hands during operation. Press the "Start" button to begin operation. The default session time is set for 30 minutes, but you can adjust the duration using the "Time +" or "Time-" buttons as needed.
4. You can adjust the energy intensity during use to enhance comfort, using the"Energy +" and "Energy-" buttons.
5. The operational functions of this device can be controlled both from the screen on the machine and via the provided remote control, which is dedicated to this specific unit. After use, turn off the power switch and disinfect the device with 75% alcohol to prepare for next use.

6. Instructions for using the belt:

First, connect the Belt to the device and secure it to the appropriate part of your body, ensuring it fits snugly against the body (it can be worn over thin clothing). Once connected to the device, turn it on to begin use.

Magnetic vibration function

Press the "Magnetic Vibration" button on the screento enter the "Low" magnetic vibration message mode.Press the "Magnetic Vibra-tion" button again to switch to the "Medium" or "High" massage settings.

Magnetic heat function

Press the "Magnetic Heat" button on the screen to enterthe "Low" magnetic heat mode. Press the "Magnetic Heat" button again to switch tothe "Medium" or "High" magnetic heat settings.